

**Instruction Manual
BS-8042**

Bath Electronic Scale

Bath Electronic Scale BS-8042

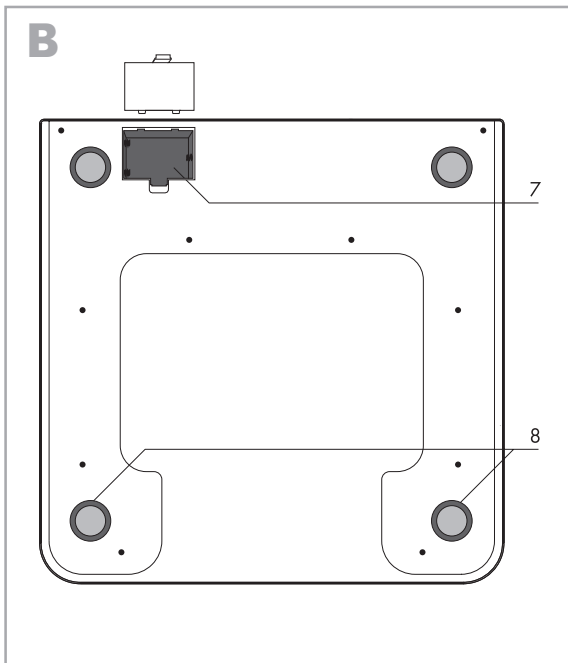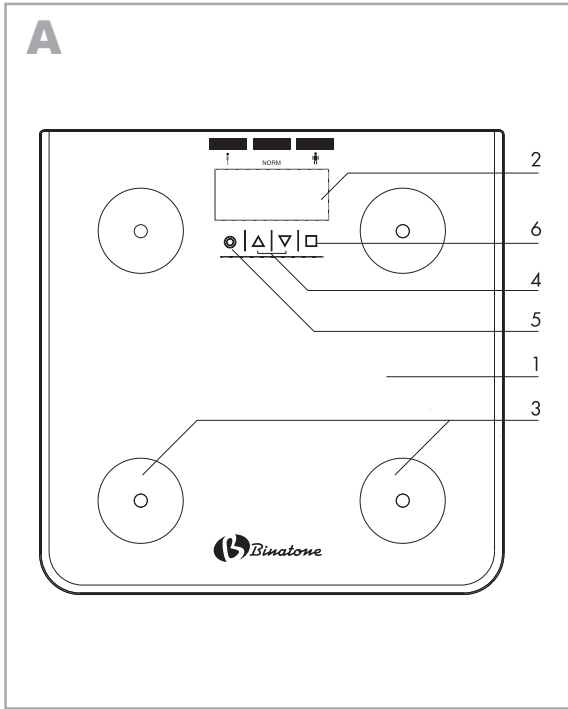ENGLISH4

ENGLISH

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Bath Electronic Scale BS-8042

DESCRIPTION OF PARTS OF THE APPLIANCE

Picture A:

1. Tray
2. LCD display
3. Electrodes
4. Up/Down buttons
5. On/Set Button
6. Confirmation button

Picture B:

7. Battery compartment
8. Legs of the appliance

BEFORE THE FIRST USE

Before the appliance is used for the first time, remove all packing. Insert the batteries into the battery compartment **(7)** on the platform **(1)**.

Be careful! Do not use the appliance if you have cardio stimulators or other electronic medical appliances inside! Make measurements at the same time of the day as the appliance is sensible for body water content. Always step on with bare feet in the middle of the platform with the electrodes.

OPERATION

Weighing

1. Place the scale on a flat solid horizontal surface. Do not put the scale on carpets or other soft coatings, since this will affect the accuracy of measurements.
2. To switch on the scale, press On/Set button **(5)** on the platform **(1)** and wait till "0.0" appears on the display. Choose the measurement — kg, lb, st by pressing Up/Down buttons **(4)**. Evenly distribute the weight, and stand still. Keep waiting for several seconds, until the appliance counts the result, which will first blinks two times on the display **(2)** and then fix up.
3. When measurements are over, get off from the platform **(1)**. The scale will automatically switch off in 8 seconds.

Note 1: *if the weight measured exceeds the permissible load (180 kg), the display **(2)** will show «OL» indicating that the permissible load has been exceeded.*

Note 2: if the scale display (2) will show «Lo», this means that the batteries are exhausted and need replacement. Open the compartment cover (7) and carefully remove the batteries. Observing the polarity, insert new batteries in the appliance.

Note 3: if after the weighing, the display (2) shows «Err», this means that an error occurred during the weighing process. Repeat the weighing procedure again. Measurement errors can occur when fat percentage is under 3% or exceeds 50%.

Personal data input

Functional buttons (4-6) to input height, weight, sex and other personal data are on the platform (1).

Allowable range of parameters:

- ▶ memory: 1–10 Users,
- ▶ height: 100–220 cm,
- ▶ age: 10–80 years,
- ▶ weight: 2–180 kg,

User's number

1. Press On/Set button (5). By default, the following data will appear on the display (2) sex, height, age. At the same time, the User's number will be blinking.
2. Use the button (4) to increase/ decrease User's number.

Note: at each one-time pressing of the button (4), readings will be increased/ decreased by one.

Sex

1. Press confirmation button (6) to confirm the User's number and to enter the User's sex input mode.
2. Use the button (4) to change the User's sex (icon of female/ male).

Height

1. Press confirmation button (6) to confirm the User's sex and to enter the User's height input mode.
2. Use the button (4) to increase/ decrease the User's height (from 100 to 220 cm).

Age

1. Press confirmation button (6) to confirm the User's height and to enter the User's age input mode.
2. Use the buttons (4) to increase/ decrease the User's age (from 10 to 80 years).

Note: if you do not introduce any datum within 8 seconds, the scale will switch off automatically.

Body fat/water content measurement

1. Switch on the scale by pressing On/Set button (5). The LCD (2) will show the data set by the last weighted User.

2. By pressing the up/ down buttons **(4)**, select your User's number with the personal data you introduced.
3. Before weighing, wait till the weight reading indication «**0.0**» appears on the display **(2)** .
4. Now you can stand up on the platform **(1)** and weigh yourself. Keep waiting for a few seconds to allow the appliance to calculate the result.
5. Once finished the weight, body fat percentage, water percentage, muscle percentage, bone weight, kcal and Body Mass Index will be displayed alternatively for 3 times before the scales automatically turn off.
6. One of 3 indicators will light on showing the fat percentage (blue – underweight, green – healthy, red – overweight).
7. The scale switches off automatically once the measurement is over.

Note 1: if you step on the electrodes **(3)** wearing shoes or stocks, the scale will not be able to calculate the percentage content of water/ fat, and in this case the display **(2)** will show “---”.

Note 2: body water and fat readings and recommended calorie intake will appear 3 times one after another.

Note 3: if the body fat content is lower than 3.0% and exceeds 50.0%, “Err” will appear on the display **(2)**.

Table of the fat/water ratio percentage content measurement results

Having received the weighing results and measured the body fat/water content, check the state of your organism against the values indicated in the Table below:

Age	Fat content (female)	Fat content (male)	Status
<=30	3.0–16.4%	3.0–13.6%	Underweight
	16.5–22%	13.7–18.0%	Healthy
	22.1–30.4%	18.1–26.4%	Overweight
	30.5–50.0%	26.5–50.0%	Obese
>30	3.0–17.6%	3.0–14.4%	Underweight
	17.7–23.0%	14.5–19.0%	Healthy
	23.1–32.0%	19.1–28.0%	Overweight
	32.1–50.0%	28.1–50.0%	Obese

Water/Muscle Ratio Chart

Age	Water content (female)	Muscle (female)	Water content (male)	Muscle (male)	Status
<=30	53%	44%	56%	46%	Underweight
	53-58%		56-61%		Healthy
	58%		61%		Overweight
>30	51%	44%	54%	46%	Underweight
	51-56%		54-59%		Healthy
	56%		59%		Overweight

Bone Weight Chart

Body weight/bone weight (female)	Body weight/bone weight (male)
45 kg/3 kg	60 kg/4.5 kg
45–60 kg/4.2 kg	60–75 kg/6.0 kg
60 kg/6.5 kg	75 kg/7.5 kg

Kcal Chart

Age	Female		Male	
	weight	kcal	weight	kcal
6–17	50	1265	60	1620
18–29	55	1298	65	1560
30–49	60	1302	70	1561
50–69	60	1242	70	1505

BMI Chart

BMI	LOW	HEALTHY	HIGH	OVERWEIGHT
		18.5	18.5–27	30–39

CARE AND CLEANING

1. Wipe the appliance with a damp lint free cloth, wiping it dry, so that there are no water marks on the tray **(1)**. Do prevent water getting inside the appliance. Do not use hard abrasive materials and substances for cleaning.
2. Always keep the appliance in a horizontal position in a dry, cool place. Do not store the appliance in places with high humidity or temperature.

SAFETY MEASURES

This appliance is powered by batteries, therefore during operation it is necessary to comply with the following safety measures:

- ▶ Before using the appliance with batteries, study directions for their use.
- ▶ Only use the batteries as recommended by the manufacturer
- ▶ It is not recommended to use batteries and accumulators in the appliance at the same time, since it will cause quick discharge of the battery.
- ▶ It is strictly prohibited to re-charge batteries as it can cause an explosion.
- ▶ To charge re-chargeable batteries use only chargers recommended by the manufacturer.
- ▶ Before turning the appliance on make sure the batteries are set correctly and are in good working conditions.

- ▶ Do not use batteries with visible marks of oxidation.
- ▶ If the corrosion liquid comes in contact with the skin or the eyes, immediately wash them with plenty of water, and if any is ingested – drink a lot of water, cause vomiting and go to the doctor.
- ▶ If the appliance is not going to be used for a long period of time, it is recommended to take the batteries out of the appliance.
- ▶ Do not use the appliance where the air may contain volatile flammable substances.
- ▶ In cases where liquid has got inside the appliance, take out the batteries. See the attached instruction to become acquainted with their further use.
- ▶ Do not use the scale near electromagnetic radiation sources (computers, radio telephones etc.).
- ▶ Do not keep the scale near heating appliances, heating systems, or in the direct sun.

SPECIAL SAFETY MEASURES

- ▶ Use the appliance only for purposes stipulated by this instruction book.
- ▶ Do not try to repair the appliance independently, go to an authorized service center.
- ▶ Only perform the weighing of children and people with limited physical abilities under strict supervision.
- ▶ Only use the scale on a hard horizontal surface.
- ▶ Do not stand on the tray with wet feet since it can be slippery. Always make sure that the scale is dry.
- ▶ While weighing, always stand in the central part of the tray in order not to overturn the scale or fall off.
- ▶ Try to evenly distribute the load on the tray of the scale.
- ▶ During weighing do not lean over other objects.
- ▶ Whilst standing on the scale, do not make quick movements, jumps or kicks, as this will cause incorrect operation of the scale or breakage to the glass and plastic.
- ▶ Only weigh within the limits, indicated on the box of the appliance.
- ▶ Weighting is recommended to perform at the same time.
- ▶ Results of measuring of fat/water content in the body are not a medical purpose and are only for the reference.
- ▶ Do not insert foreign objects in the open parts of the appliance.
- ▶ Prevent the appliance from falling or other objects falling onto the appliance. As this can damage the appliance and affect the accuracy of measurement during weighing.
- ▶ Prevent water getting inside the appliance.
- ▶ Do not use the scale or store them in places with high humidity or where water could be spilled on them.
- ▶ Store the appliance in the original package, in a horizontal position.
- ▶ Do not put foreign objects on the scale.
- ▶ When cleaning the scale, do not use abrasive materials and cleansers as these can damage the body and parts of the appliance.

AFTER-SALES SERVICE

If your appliance does not operate or does not operate properly, contact an authorized service center for consultancy or repair.

The warranty does not cover:

- product packing, consumables and accessories (knives, tips, etc.);
- mechanical damages and natural wear of the product;
- failures due to overloading, improper or careless maintenance, penetration of liquids, dust, insects and other foreign matters inside the product, exposure of plastic and other not thermal-resistant components to high temperatures, force-majeure circumstances (accidents, fire, flooding, power grid failure, lightning stroke, etc.);
- damages caused by:
 - ▶ Non-fulfillment of the above-mentioned precautions and operating manuals.
 - ▶ Application in commercial purposes (application of the product beyond personal daily living needs).
 - ▶ Unauthorized repair or replacement of the components of the appliance, performed by an unauthorized person or service center.

Also the warranty doesn't cover replaceable batteries.

ENGLISH

DISPOSAL

This appliance should be disposed of in accordance with waste disposal regulations in your country and with the least damage to the environment.

Only ecologically safe materials are used in the packaging of our products and the cardboard and paper can be disposed of with normal paper waste.

Manufacturing date is specified on the appliance.

RU XX YY ZZZ code reads as:

XX – production year

YY – production month

ZZZ – batch number

Major Importer of Binatone TM in Russia: JSC Interworld Products. The address: 2-nd Khoroshevsky proezd, 7, b.1, Moscow, Russia, 123007.

Manufacturer: Binatone Industries Ltd, Great Britain, U.K., 1 Highbridge Wharf, Greenwich, London, SE 10 9PS.

Made in PRC for Binatone Industries Ltd, Great Britain.

Exploitation period: 3 years after the day of purchase.

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